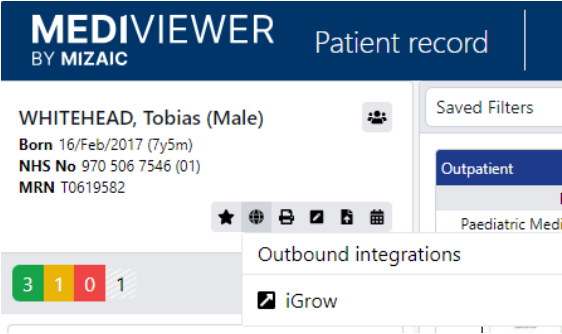
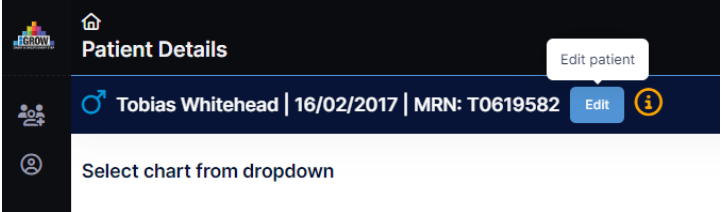
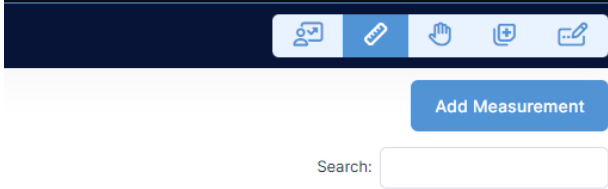
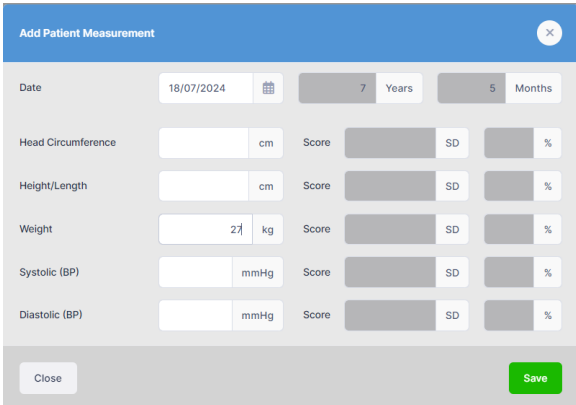
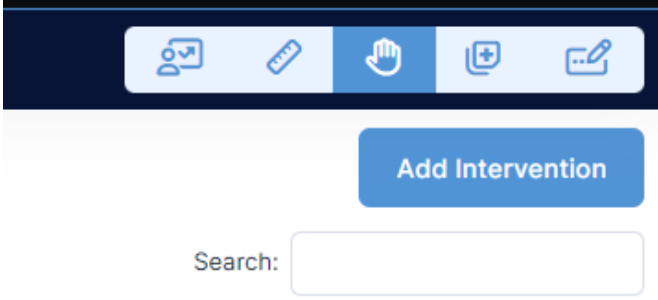
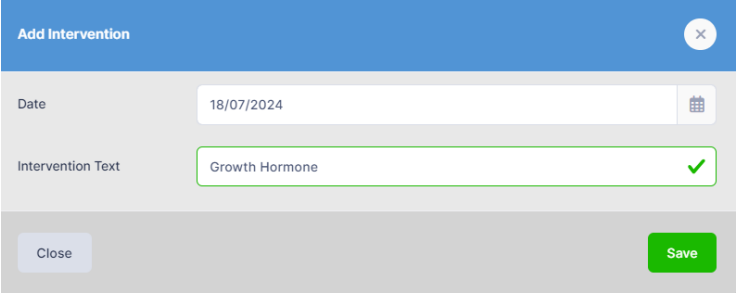
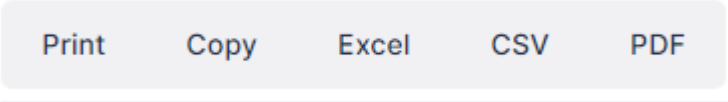
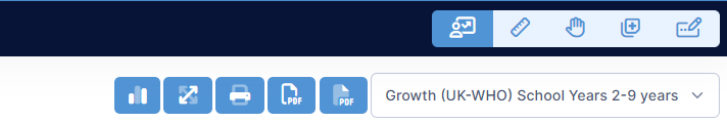

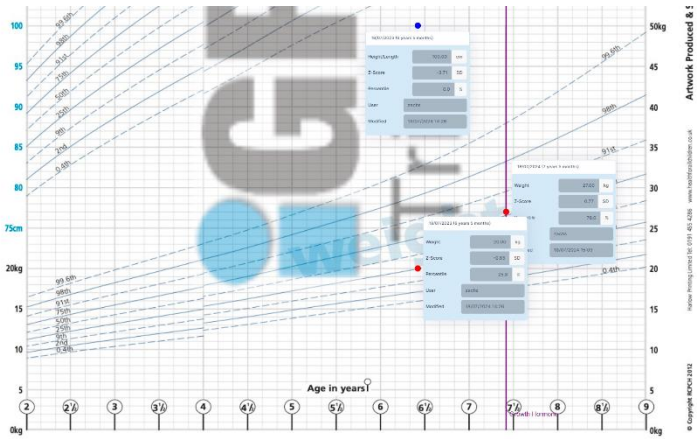


iGrow	
<p>You can access iGrow by opening a patient in Mediviewer and selecting the 'Outbound Integration Menu' from their banner and then selecting 'iGrow' from the integrations list</p>	
<p>This then opens their iGrow record. When opening the record, it is best practice to first check the patient information and if not already added and relevant, input the estimated due date and/or syndromes as this will have an effect on which charts are available and the way data is plotted (e.g. adjusted age for prematurity, Down's syndrome growth chart).</p>	
<p>To add a measurement, click on the ruler/measurement icon on the top right of the record and then click 'Add Measurement'</p>	
<p>Input the data as appropriate and click Save</p>	

iGrow	
<p>To add an Intervention, click the Hand icon/Intervention tab on the top right of the record and then click 'Add Intervention'</p>	
<p>Input the Intervention in the free text box, check the date is correct and click Save</p>	
<p>If needed you can Print, Copy or Export Measurements and Interventions to Excel CSV or PDF Files by selecting the relevant button in either tab.</p>	
<p>Select Chart View by selecting the Charts Tab on the top right of the record.</p> <p>You can then use the dropdown box to select the chart that you wish to look at.</p>	

iGrow	
<p>With a Chart selected you can use the following options (from left to right):</p> <ol style="list-style-type: none"> 1. Chart Measurements (Quick glance view of measurements on current chart) 2. View Full Screen (Scroll and click on data points for more information) 3. Print current chart 4. Export current chart to PDF 5. Export All charts to PDF 	
<p>By clicking View Full Screen on a chart we can scroll through the Chart, click on coloured Datapoints to show more information about the measurements and view Interventions (shown as a labelled purple line through the chart)</p>	
<p>To exit a full screen Chart, scroll to the top and click the close icon.</p>	